



HealthAtlas

Detailed User Manual Local Apple Health visualisation for macOS

English

Privacy-first local Apple Health visualisation for macOS

HealthAtlas - User Manual

Welcome

HealthAtlas reads an Apple Health export that you choose locally on your Mac. You then decide exactly which recognised data types appear in the dashboard. The app does not diagnose conditions or recommend treatment.

- Every normal launch starts without imported health data.
- Imported data stays in memory for the current session only. It is not uploaded or sent to an account.
- Language, theme, card count and data-type selection are stored separately for the Dev, Beta and Final app variants.
- While no personal data is loaded, the start screen also offers the matching GitHub manual and optional AI help. It copies only a general question with the public manual link to the clipboard; health values and local data stay in the app.
- Requirement: macOS 26 or later.


Important

HealthAtlas describes values and trends. It is not a substitute for professional medical advice, examination or diagnosis.

Quick start in 5 steps

The repository includes a fully synthetic demo export for a safe first test. It contains no personal health information and fictional values for every currently supported, non-deprecated Apple Health export type.

- Open HealthAtlas. The empty overview presents a central import button.
- Click "Import ZIP or Export.xml ...".
- Choose Demo/AppleHealthDemo/Export.xml or your own Apple Health export.
- Open Sources after import and turn the data types you want on or off.
- Explore the selected data in Overview, Trends and Insights.

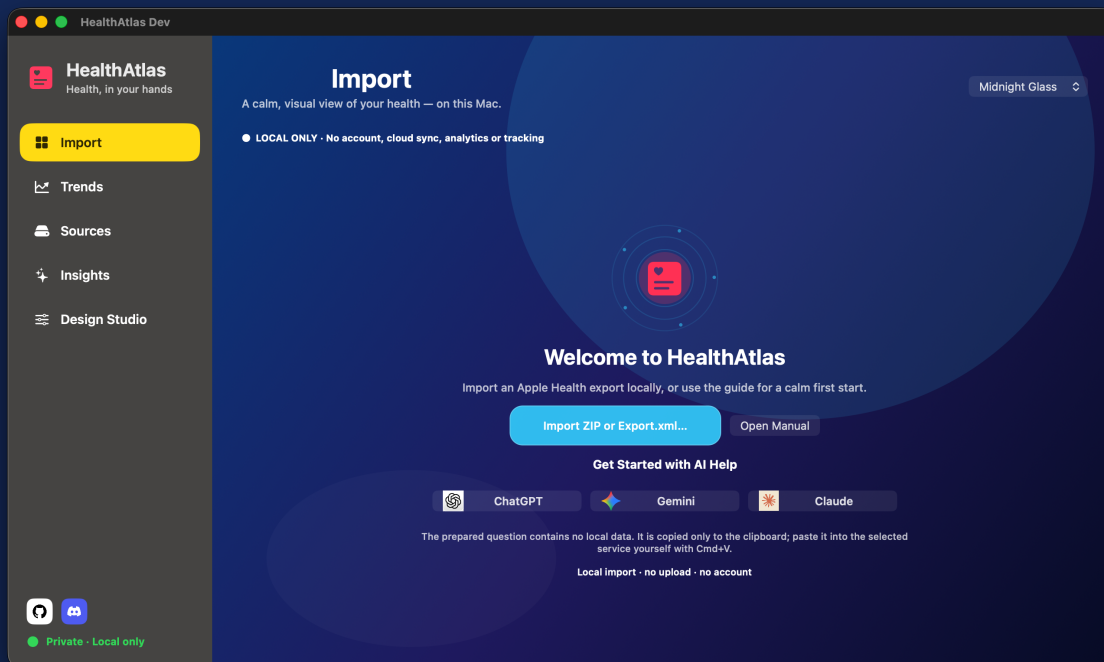

Empty start screen. The shown interface contains no health information.

Importing an Apple Health export

HealthAtlas accepts exactly two local formats: a direct Export.xml file or an Apple Health ZIP archive containing Export.xml. You do not need to unpack the ZIP first; HealthAtlas reads Export.xml directly from the archive.

- Files must be between 1 byte and 100 MB.
- Only one file can be selected at a time.
- The clinical companion file in an Apple Health export is intentionally not imported.
- There is no direct HealthKit connection and no cloud import.
- For a wrong ZIP, unreadable XML or an oversized file, HealthAtlas explains the issue and imports nothing.

Export on iPhone

In Apple Health, open Summary, tap your picture or initials, then select "Export All Health Data". Apple can change exact interface labels. Source: Apple Support, support.apple.com/en-in/guide/iphone/iph5ede58c3d/ios

Navigation and status


The frosted sidebar is visible by default in every theme. Before the first local import, its first entry is called "Import"; after a successful import it becomes "Overview". It shows "Private - Local only" at the bottom as a permanent privacy reminder.

- Overview: cards for selected data types.
- Trends: an interactive timeline for one selected type.
- Sources: selection of all data types recognised in the import.
- Insights: a descriptive summary for one data type.
- Design Studio: language and appearance.
- The native menu bar provides import, the local PDF report, Design Studio and window controls. Use View to show or hide the sidebar.
- Interactive cards and charts are available as controls for macOS assistive technologies.
- The GitHub and Discord icons above the privacy status open the project page and community in the default browser.

Overview

After a successful import, Overview shows only the types enabled in Sources. Every card has a type-specific accent colour and graphic treatment. The displayed numbers use the most recent available daily value.

- Your period at a glance only describes how many enabled types contain values across the latest seven locally recorded days — never a goal, rating or diagnosis.
- Choose 4, 8 or 12 to control cards per page. Compact, Standard and Focus change card density locally. Drag and drop cards to reorder them; each layout saves its own order and can be reset.
- Previous and next controls appear when more selected types exist than fit on one page.
- Click a card to open a true full-screen focus view with a trend, period comparison and yearly calendar; Open Trends leads to the detailed view. Exit full screen returns to the app. Before a local PDF report, choose period, data types and report theme; it writes only to a location you choose.
- Below the cards, a shared multi-metric timeline and daily rings show recent local values from several types. Rings are neither goals nor ratings. At narrower window sizes, legends remain within their card.
- Steps, energy, distance and flights climbed are shown as daily sums. Other numeric types are shown as daily averages.
- Name, value, unit and date come from the selected export and follow the app language formatting.

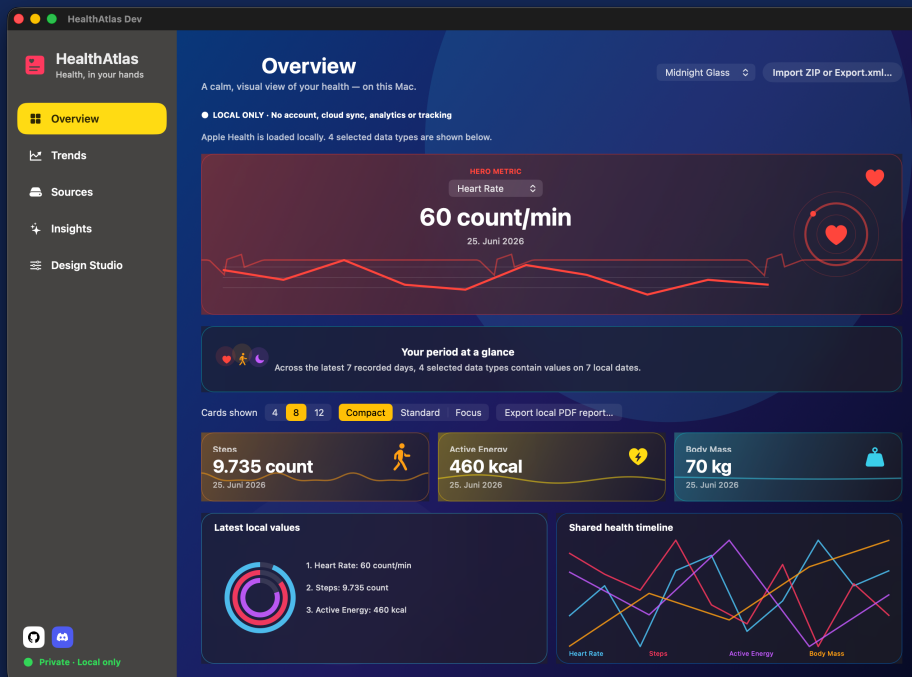


Overview with synthetic demo data. Colours aid recognition; they do not assess health.

Sources

Sources is the central selection area. It lists every recognised data type with an on/off switch, number of samples and a summary value. The table scrolls vertically when an export contains many types.

- Show all enables every recognised data type.
- Show none disables all data types.
- The switch in each row immediately determines whether a type is available in Overview, Trends and Insights.
- Search and categories help with large exports. First choose Overview, Trends or Insights as the pin target; the star then marks a favourite only for that area. Arrows set the order for the current dashboard layout.
- For a first personal import, the first four recognised types are selected by default. A compatible earlier selection is reused when its types occur in the new import.
- Local data quality counts selection, dated types, local days and samples. It also names missing days in the current local period and types with few values. It does not rate health data.
- Replace import opens the local file picker again. Delete all local data clears the current session only after confirmation. The visible import time applies only to the current session.

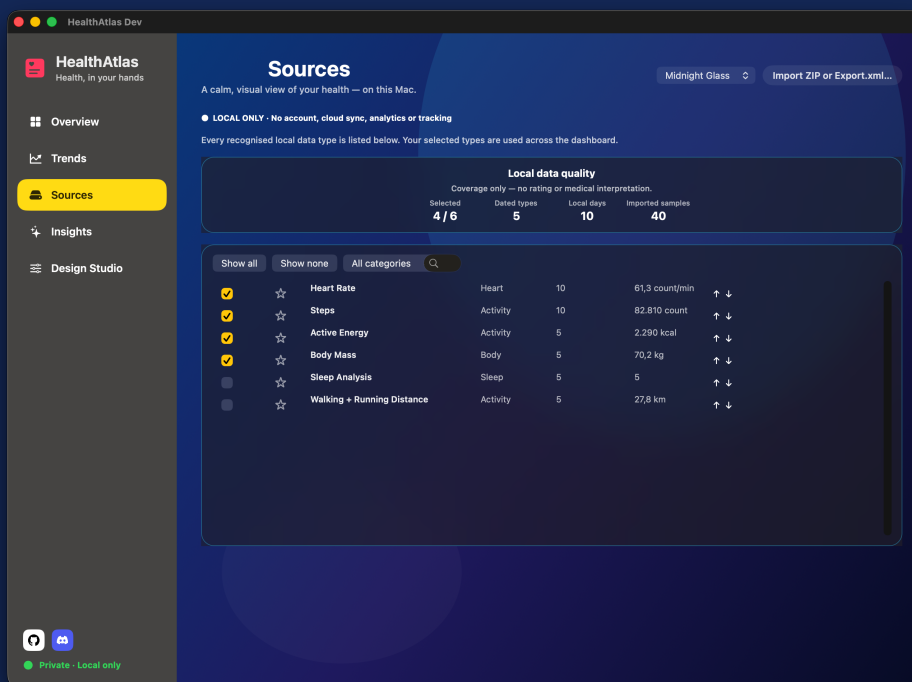

Sources with synthetic demo data. This selection controls the content of all other areas.

Trends

This section visualises one active type according to its local format: bars for steps and energy, an area for sleep, and a line for other numeric types. Select a period first, then choose one of the enabled types in the menu.

- 7D: last 7 days.
- 30D: last 30 days.
- 3M: last 90 days.
- 1Y: last 365 days.
- Hover a point to reveal a floating value card with date, formatted value and mini-trend. Clicking also gives that point a gently pulsing marker. Factual highlights can name the highest local value in the period or the number of dates with data.
- A local comparison places the current and immediately preceding periods side by side. It is not a rating or diagnosis.
- The local data calendar sits in a separate panel below the chart so day cells and line animation remain clearly separated.
- If fewer than two values exist in the selected period, the chart reports that there are not enough values.




Trend with synthetic demo data. Points are interactive and reveal their value after a click.

Insights

Insights condenses one enabled data type into a local descriptive snapshot. Choose the type from the menu. The view shows the latest daily value, its date, local coverage and a recording pattern; it intentionally contains no detailed trend or data calendar.

- Local coverage only names days with and without a value in the current short period. The pattern names the most frequently recorded weekday. Neither is a rating.
- Local pattern only describes which weekday has the most locally recorded dates — not health or behaviour.
- The screen does not supply normal ranges, alerts, ratings or medical conclusions.
- If only one dated value is available, the app states this clearly.

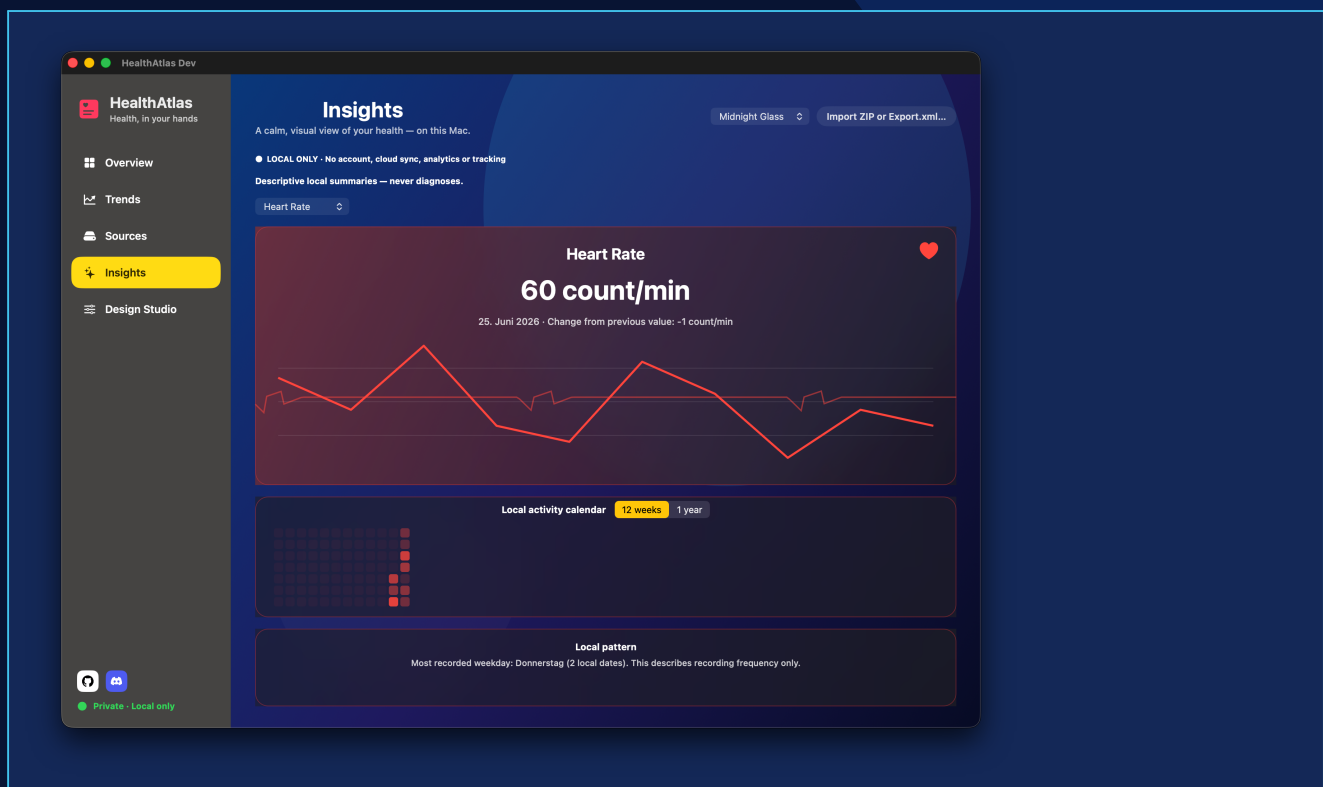


Insights with synthetic demo data. The card is a descriptive summary, not a diagnosis.

Design Studio

Use Design Studio to set the interface language and appearance. Changes take effect immediately and are stored locally for the current app variant.

- Language: Deutsch or English. Navigation, labels and known data-type names change with the selection.
- Clear Glass: a shared frosted layer with a calm cyan, blue, violet and pink glow plus subtle light points behind cards, text and controls. This keeps values readable throughout the app.
- Midnight Glass: a dark blue glass surface.
- Aurora: a teal variation.
- Warmpaper: a warm red-violet variation.
- Cards enter in a staggered sequence, react subtly on hover, and import completion briefly shows a success shimmer. Changing a range or data type redraws the chart softly. Page and sidebar changes use gentle transitions.
- Clear Glass adds very subtle moving contours and transparent health symbols to its glow and sparks. Reduce Motion and importing substantially reduce animation.
- When macOS Reduce Motion is enabled and while importing, Clear Glass reduces or pauses motion.
- HealthAtlas opens in 16:9 and remains freely resizable afterwards. The layout adapts to the selected window size.

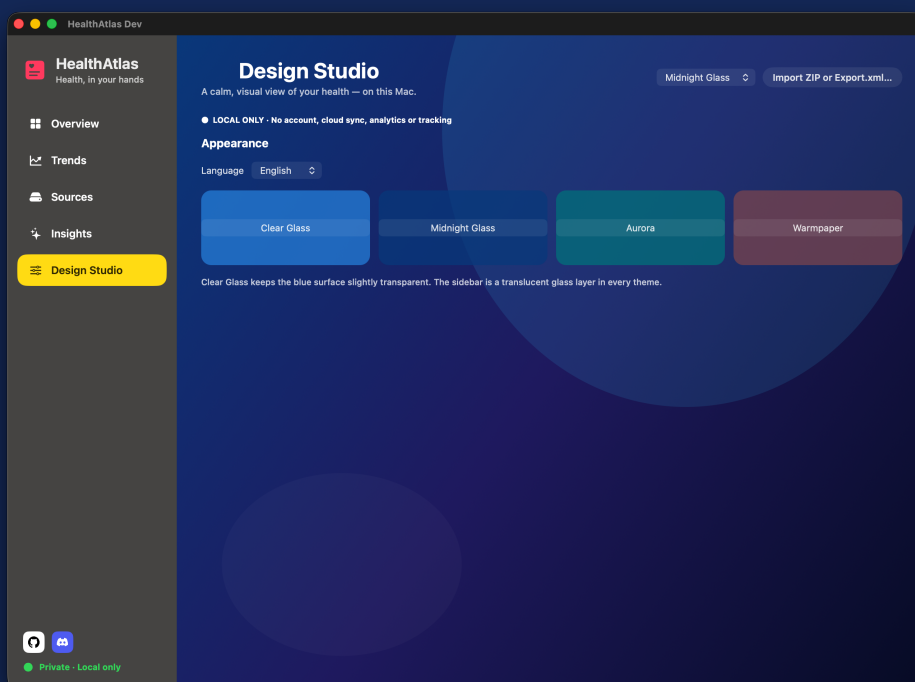

Design Studio. Themes change appearance only, never the health data.

Privacy and limits

HealthAtlas is designed as a local visualisation. It has no account, analytics, advertising, tracking or hidden upload. The app does not send imported health values to HealthAtlas, GitHub, Discord or another service.

- Data remains in memory while the app is open and is not loaded again at the next normal launch.
- GitHub and Discord are opened only by an explicit click on the lower sidebar icons. They receive no imported values.
- ChatGPT, Gemini or Claude open only after an explicit click on the empty start screen. Before that, HealthAtlas only copies a fixed general introduction question with the public, language-specific manual link to the clipboard. Only Cmd+V lets the person decide whether to paste the question into that service.
- The bundled demo is synthetic and allows feature testing without personal data.
- HealthAtlas is not medical software and not emergency assistance.

Symptoms or uncertainty

Please contact qualified healthcare professionals. A visual change in a chart is not a medical statement.

Troubleshooting and Gatekeeper

HealthAtlas is ad-hoc signed because the project has no Apple Developer account. macOS Gatekeeper may therefore warn when opening a Dev, Beta or Final build for the first time.

- In Finder, Control-click the app and choose "Open". Confirm Open in the following dialog.
- If necessary, go to System Settings > Privacy & Security and choose "Open Anyway" for that exact HealthAtlas build.
- Do not disable Gatekeeper system-wide. Open only your own build or one from the official HealthAtlas project.
- For import errors, confirm that the ZIP contains Export.xml, the file is not empty, is below 100 MB and is stored locally.
- If no values appear, enable at least one type in Sources. Trends need at least two daily values in the chosen period.

Build variants

Dev, Beta and Final use separate local preferences. A theme or selection in Dev does not alter the settings of a Beta or Final app.

Complete control reference

This table summarises every current area, control and outcome.

Area	Control	Outcome
Header	Theme menu	Changes appearance immediately.
Header	Import	Opens ZIP or Export.xml picker.
Menu bar	File / View / Window	Imports, exports the local PDF report, opens Design Studio, shows or hides the sidebar, and controls the window.
Overview	Period	Describes only local recording days.
Overview	4 / 8 / 12 · density	Sets card count and local card density.
Overview	Arrows	Moves through additional selected cards.
Overview	Card / PDF report	Opens the focus view or saves a local PDF report.
Overview	Timeline / rings	Shows several local types; no goal or rating.
Trends	7D / 30D / 3M / 1Y	Limits displayed days.
Trends	Data type menu	Changes the displayed timeline.
Trends	Data point / hover	Highlights a point or shows date, value and mini-trend.
Sources	Row switch	Enables or removes one type.
Sources	Search / category / star / arrows	Filters, favourites and orders types locally.
Sources	Show all / none · data quality	Enables types or counts local coverage only.
Insights	Data type menu	Changes the descriptive snapshot.
Design Studio	Language	Switches Deutsch and English.
Design Studio	Theme cards	Selects Clear Glass, Midnight Glass, Aurora or Warmpaper.
Lower sidebar	GitHub / Discord	Opens the external project or community link in the browser.